



**brite brush**  
sonic on the go




**SONIC ON THE GO**



# SONIC ON THE GO

## Easy Brushing

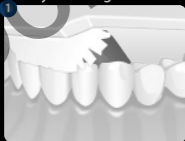

Soft bristle design for easy brushing, brushes hard to reach places and dramatically reduces plaque.

## Inserting Batteries



Turn bottom battery cover counter clockwise. Insert (1) AAA battery. After inserting battery replace battery cover and power on to start brushing.

## Replacing Brush Head



### Removing Old Brush

Hold your brite brush sonic on the go toothbrush firmly. Rotate brush head counter clockwise and pull brush head upwards.

### New brush head


Insert new brush head on toothbrush base and turn counter clockwise until brush head locks into place.

## Powering ON your brush



Press power on button as showed in picture. Make sure to rinse your brite brush sonic on the go thoroughly after each use. Cover toothbrush each time with the sanitization cap for maximum sanitization.

## Sonic on the Go

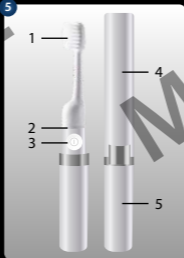

1. brite brush sonic on the go promotes highly efficient sonic technology for maximum sanitization.
2. Replaceable brush heads
3. Power On/Off button.
4. Sanitization Cap
5. Battery Cover

## Medical Warnings

1. Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
2. Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
3. If you have a pacemaker or other implanted devices contact your physician or the device manufacturer prior to use.
4. Consult your physician prior to using the Brite Brush sonic on the go® if you have medical concerns.

